

VascuPen 0.5 – 2 mm Tips

Indications:

- Skin Vascular Lesions

Pre-Treatment:

- Clean treatment area; remove all make-up or moisturizing cream.
- Always perform a skin test on the intended treatment area using the lowest recommended fluence before the first treatment session. Assess skin condition after 15 – 30 min.
- All personnel within the laser operating area must wear appropriate eye protection.

Treatment:

- Position the applicator perpendicular to the skin over vascular lesion.
- Pulses should be gapped at least 0.5 – 1 cm apart. Do not perform a second pass.
- Endpoint: change in vessel coloration (reddening, darkening, or bleaching)
- Number of treatments: 1 – 3 sessions.
- Treatment intervals: 6 – 8 weeks.



Recommended Treatment Parameters – Skin Vascular Lesions:

Fitzpatrick Skin Type	Treatment Area	Spot Size (mm)	Fluence (J/cm ²)	Pulse Width (ms)	Frequency (Hz)
I – IV	Very small vessels	0.7	18	19	1
	Vessels	1	18	23	1
	Vessels	2	20 - 22	51-63	1

Post-Treatment:

- Avoid sun exposure after and in between treatments, use UVA/UVB sunscreen (SPF>30).
- If treatment area shows erythema, apply Trolamine, Cicaplast or equivalent skin calming products.

VascuPen 1-2 mm Tip

Indications:


- Warts

Pre-Treatment:

- Clean treatment area; remove all make-up or moisturizing cream.
- Always perform a skin test on the intended treatment area using the lowest recommended fluence before the first treatment session. Assess skin condition after 15 – 30 min.
- All personnel within the laser operating area must wear appropriate eye protection.

Treatment:

- Position the applicator perpendicular to the lesion.
- Use a stamping technique over the treated area with no overlap.
- Endpoint: Lesion darkening or detachment.
- Number of treatments: 1 – 3 sessions.
- Treatment intervals: 6 – 8 weeks.

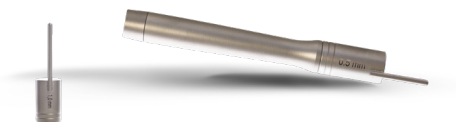

Recommended Treatment Parameters – Warts (1 – 2 mm tips):

Fitzpatrick Skin Type	Spot Size (mm)	Fluence (J/cm ²)	Pulse Width (ms)	Frequency (Hz)
I – IV	1	16-20	13-23	1

Post-Treatment:

- Avoid sun exposure after and in between treatments, use UVA/UVB sunscreen (SPF>30)
- If treatment area shows erythema, apply Trolamine, Cicaplast or equivalent skin calming products.